


# Benefits

## Flexible image acquisition

- Integration of various **flat panel and CR systems** by different manufacturers
- Option to **connect up to 3 flat panels** (bucky, wall stand and mobile) to one system
- The **configurable generator interface** enables the user to control X-ray generators or X-ray systems by different manufacturers, delivering the generator settings directly from the software
- Option for the **parallel operation of a flat panel and a CR system** included in the standard package. The user has the choice to take the next image with either the flat panel or the integrated CR system. This flexibility also provides an **excellent emergency concept** in case of a defect flat panel.

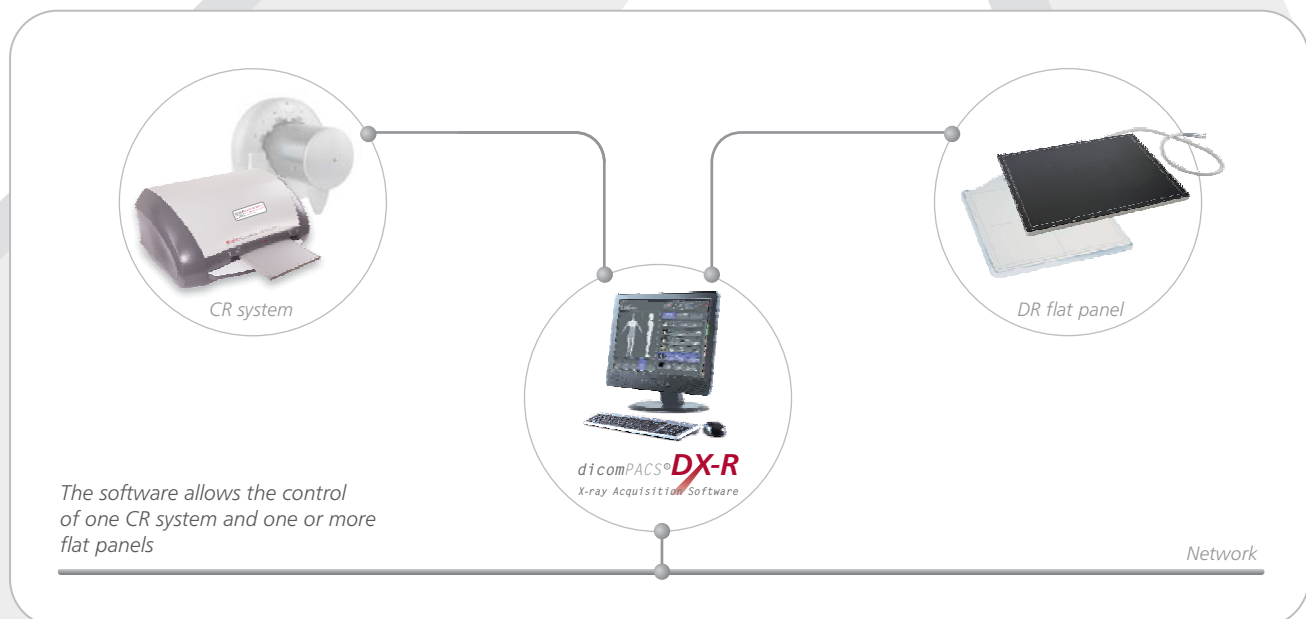

- **AEC** (Automatic Exposure Control) and **ARP** (Anatomical Programmed Radiography) allow the user to **automatically adjust all X-ray options** for each projection with an option to subsequently edit the image manually
- Integration of **dose area product meters** (DAP) – the readings are saved directly to the relevant image
- Electronic X-ray log