



After the patient has been called up from the electronic card index, the preoperative images can be displayed immediately. If individual images are to be used in the operation report to be created later, they can be identified by pressing a button and so saved in the operation programme.

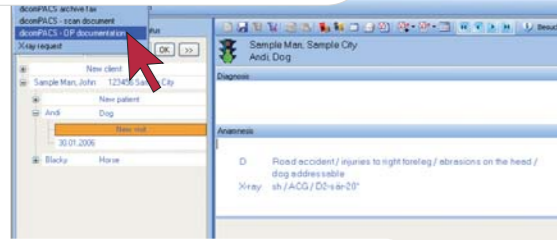

1. Patient called up from management system

A click on the Record button starts the image recording.

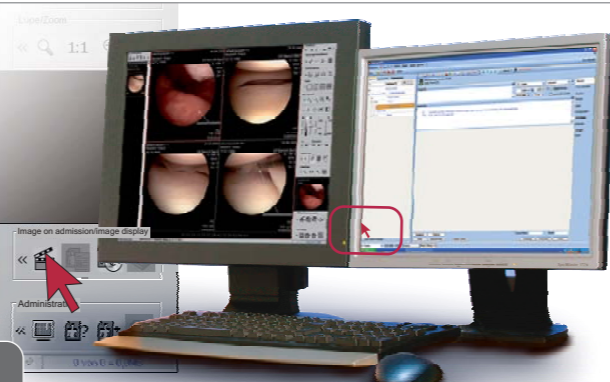


2. Start recording

By means of a foot switch or the endoscope trigger the individual images or video clips are recorded. It is possible to switch between image sources during the operation, e.g. to the C-arm, operation camera etc.



3. Performing the operation

The operator skims over the recorded images at the PC. Those intended for the operation report can then be selected and text and arrows can be added in. Images can also be imported from digital cameras.

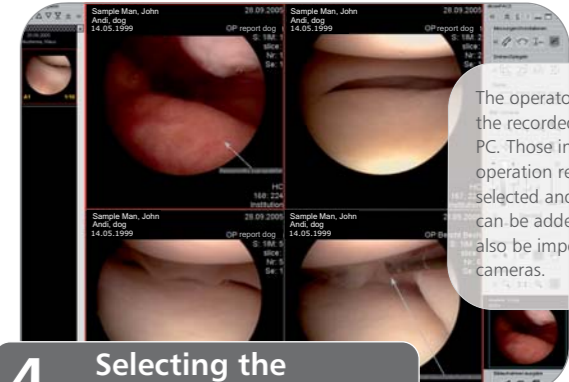


4. Selecting the operation images

Step by step

Simple user friendly workflow

5. Choosing a document template

It is now possible to select a suitable operation report template from a number of alternatives, according to the nature of the operation.

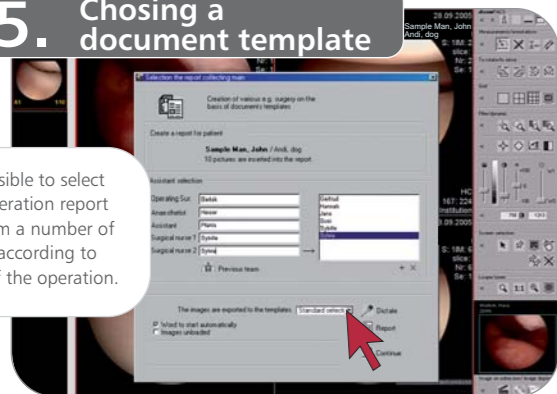


6. Digital dictation

If needed, digital dictation can be recorded during the operation. The operation report is then saved together with the images and any spoken dictation.




7. Operation report list

All prepared operation reports including the dictation are available under MS Word in a job list.

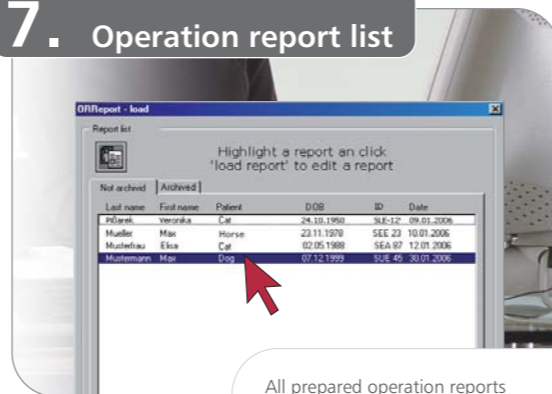

8. Creating the operation report

After clicking on the report to be completed, it is opened with its integrated images. All data about the patient, on operation images, as well as information about the operation - such as assistant or anaesthetist are automatically added in. If a dictation was made, this can be listened to and written directly into the report. This function is available at every workstation (with MS Word capability). A secretary can therefore also complete the report and print it out.

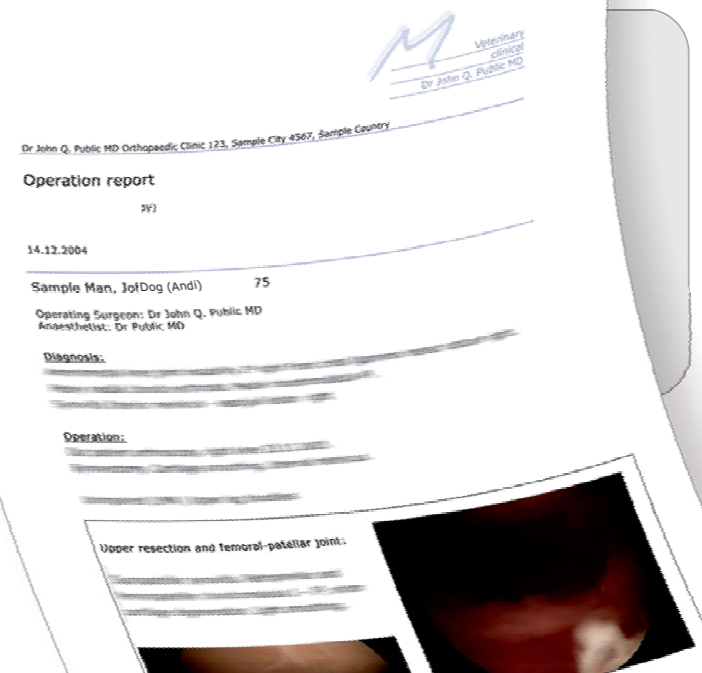

9. Printing and archiving

The operation report is now printed and automatically archived together with the examination images as a unit of image and diagnosis. It may now be called up from the patient register.

