



# Product·Information

## X-ray table CombiVet S, variant 1,2 and 3

Stand for X-ray unit with integrated holder for cassette and grid

### Brief description

- The unit support pole moves along the table, simultaneously with the built in grid ready cassette holder. The X-ray unit is always perfectly aligned with the X-ray cassette. Only minimum repositioning of the patient is required.
- Table frame: anodised aluminium
- Table top: radiolucent white synthetic
- Easy to clean and disinfectant proof
- Optional: Adjustable positioning aide for dysplasia coxae
- Dimensions:  
 Table top: 155.5 x 67.5 cm  
 Working height: 73 cm  
 Total height: 182 cm




(Without X-ray generator)




# X-ray table „CombiVet S“

X-ray unit suspension module available in three variations (Type 1 to 3)

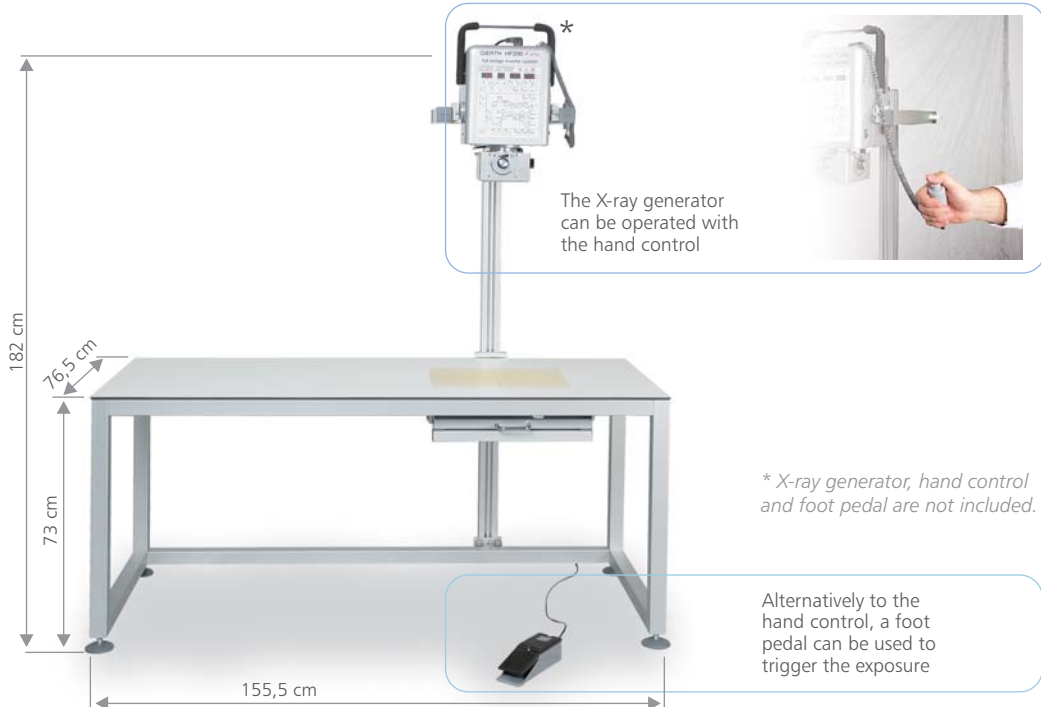


## Suspension Type I

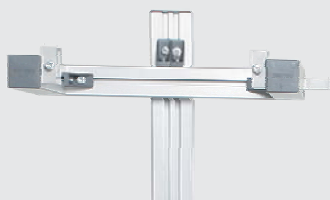

**Important: Please indicate your X-ray unit type for the required customization of the unit suspension size.**

CombiVet S **Type I** is fitted with our standard unit suspension. This suspension is not height-adjustable; the film-focus-distance is fixed.

## Suspension Type II

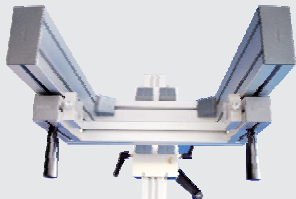

CombiVet S **Type II** is equipped with a height-adjustable unit suspension module, operated via crank case. It is also possible to slide the X-ray unit backwards and forwards in axial direction.

## Suspension Type III



In addition to height-adjustment and axial movement, a pivot joint integrated in the unit suspension of CombiVet S **Type III** allows rotation of the X-ray unit for slanted beam direction.

### Specifications subject to revision without notice

The editor strives to impart correct and up to date information. The provided specifications are based on current knowledge and are subject to revision without notice. This brochure is subject to correction. The editor assumes no responsibility for the information being up to date, correct and complete.

All furnished logos, pictures and graphics are property of the particular company and subject to copyright of the licensor. Use, dissemination, distribution or copying of the pictures, logos or text compiled or processed by the editor is subject to our written consent. All rights reserved.

Stemp of distribution partner